

Frozen Eyes Machine

User manual



✓ Please read this manual carefully before using this instrument and follow the instructions in the manual strictly.

Content

Functions.....	3
Working Principle.....	3
Product Advantages.....	3
Technical Parameter.....	4
Spare Parts.....	4
Installation Instructions.....	5
Treatment Steps.....	5
Taboo Crowds.....	7
Suggested Parameters.....	8
After Treatment Care.....	8
Machine Maintenance.....	8
After Sale Specifications.....	9

Functions

1. Ease wrinkles removal..
2. Eye bags, eye circles improving.
3. Essence introduction into skin (better absorbed).
4. Firming & tightening.
5. Eliminate eye edema.

Working Principle

Radio frequency

Radio frequency wave penetrate the barrier of epidermal basal melanocytes, shrinking the collagen fiber by heating the dermal collagen fiber to 55-65°C to tighten loose skin, stimulate collagen regeneration and restore skin elasticity gradually.

RF Technology has a strong product introduction function with special essence and skin care products to achieve a better result by using these two functions and collagen in deep skin will be growing again and the supporting structure will be reshaped so as to achieve your eyes rejuvenation.

Vibration care

20 thousands of vibration in per second to accelerate the circulation of blood, metabolism and to relieve the fatigue, tiring and swollen muscle.

Product Advantages

1. Special design for the eyes and six huge problems will be solved in one step .
2. No operation, injection and anesthesia needed which can be called non wounded and invasive treatment.
3. Safe and painless, can go back to work once the treatments ended .
4. High vibration with RF to ease your tiring eyes.


5. Mono-Polar Radio Frequency can stimulate the growth of collagen so it can bring your eyes an rejuvenation.

6. The treating area s got a roller inside to effectively massage the eyes skin and skin around.

Technical Parameter

Name of the machine	Frozen eyes
Radio frequency	2.0 MHz
Applications	General
Working mode	RF+vibration+product introduction
Treatment parts	Mainly for eyes and it can be used in face
Input power	6W-36W
Output power	0W-25W
Power	110/220V AC
Way of cooling	Wind cooling
Color	Blue to purple

Spare Parts

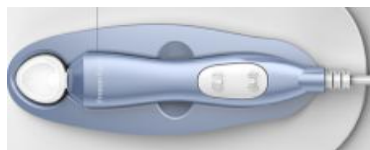

Handle



Negative clip



Power line

Installation Instructions

1. Machine assembly

- ❖ Negative clip installation: Connect the negative power line to the negative clip as the below picture shows.




- ❖ **Power line installation:** Connect the power line to the power plug back in the machine.

2. Boot


Connect the power line with machine then press the power button in the machine, and when the machine is on, the doctor can adjust the parameter for different treatment purposes.

Treatment Steps

1. Cleanse the face and take pictures before treatment.
2. Fill serum into the container in the handle;

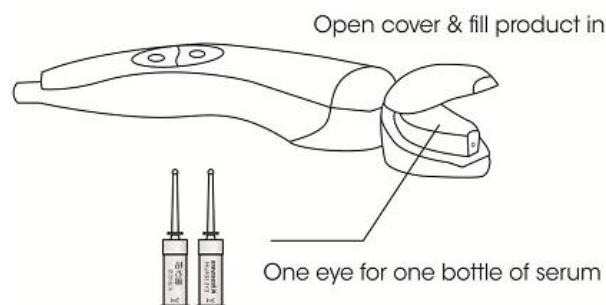

3. Negative clip usage;

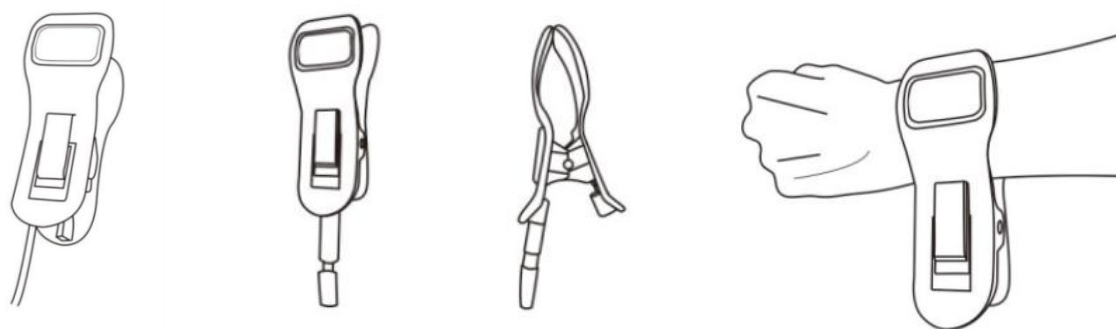

4. Treatment steps;


- 1) Circling the handle along with the eye bag and dark circle to the end of the eyebrow, repeat 3 times or more;



- 2) Move the treatment tip along with the eye bag and black circle and the fingers that are not holding the handle can move along with the tip so to lift the skin up. Repeat 5 times or more;



- 3) Circling the handle counterclockwise from the inner corner of the eye to the brow bone, then circle it counterclockwise to the end of the eyebrow and lift the skin up in the hairline area. Repeat 3 times or more;

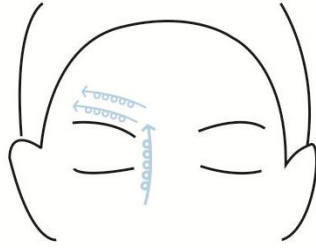

- 4) Slide from the inner corner of the eye to brow ridge and repeat the step to the end of the eyebrow so to lift the skin up and repeat 5 times or more;



- 5) Repeat the step 1-4 till 15 minutes is over and it will take about 30 minutes to finish the treatment of both eyes, Treatment course: interval time should be 7-10 days and receive the treatment 3-4 times a month.
- 6) Clean the face after treatment and take pictures for comparison.

Taboo Crowds

1. Infectious activity period and infectious disease like herpes and HIV.
2. Severely ill people(heart disease, diabetes, high blood pressure).
3. Allergies, herpes simplex and recent sun exposure.
4. Pregnancy and lactation.
5. Use other laser treatment no more than one month.

Suggested Parameters

Project	RF Energy	Vibration	Time/Minutes	Interval/days
Eye bags	2-3	1-2	20-30	7-10
Eye wrinkles	2-3	1-2	20-30	7-10
Eye fine lines	2-3	1-2	20-30	7-10
Puffy eyes	2-3	1-2	20-30	7-10
Dark circle	2-3	1-2	20-30	7-10
Face lifting	3-4	1-2	20-30	7-10
Product introduction	2-3	1-2	20-30	7-10

After Treatment Care

1. Apply the repair factor or repair eye mask right after treatment.
2. Cold compress can be applied once the treatment ends.
3. Need UV protection(stay home or wear a hat, use umbrella, sun screen).
4. Apply repairing mask or moisturizing eye mask everyday for one week after treatment.
5. Pay attention to diet, avoid spicy stimulation, seafood, beef and photosensitive food (celery, leeks, parsley, lemon), darker foods and avoid tobacco, alcohol, coffee.
6. Please do not wear makeup the day after the treatment and can wear makeup two days after the treatment.

Machine Maintenance

1. Do not place the machine in a place where exposing the machine directly to the sun and especially during summer.
2. Do not place the machine in a humid place with water vapour like sauna room or showering room or damages will do to the machine.
3. Please do not use chemical material to clean the machine or it will cause damage to the machine or machine case.
4. Do not use the machine in a consecutive time to prevent the delicate spare parts to be

damaged.

5. Please keep the machine away from children to avoid any accidents.
6. Please do not disassemble the machine case so to avoid the electrical shock and do not hit the machine or the screen or damages will do to the machine.
7. Please clean the serum on the treatment tip when finish the treatment.
8. Do not drop the handle to avoid crashing the handle.
9. Do not twist the handle cable or the handle will be broken.

After Sale Specifications

- (1) From the date of machine shipment (according to the barcode of machine factory), Party B(seller) provide Party A(buyer) one-year warranty for mainframe (excluding the handle) and lifetime maintenance service.
- (2) From the date of machine shipment (according to the barcode of machine), party B(seller) offer three months warranty for handle except the artificial factor.(artificial factor includes abnormal use of the machine, such as no water running, strong impact, fall down from high place, etc). Breaking caused by artificial factors is beyond warranty, party A(buyer) should pay the cost of breaking spare parts and two-ways freight.
 - A. Within 15 days after machine shipping from party B(seller), if handle has problem (non-artificial factor), Party B(seller) will send party A(buyer) new handle upon receipt of returned old handle (Except the nd yag or q switch laser handle, E light handle and diode laser handle); party A (buyer) bear two-way shipping cost.
 - B. During warranty (15th days to 3 months after machine shipping from party B(seller)), if handle has non-artificial factor problems, party B(seller) will offer patty A(buyer) free repairing service, party A(buyer) should bear two-way freights.
 - C. Out of warranty(more than 3 months since machine shipping from part B(seller)), if handle has any problems, Party B(seller) provides Party A(buyer) maintenance service with the cost of replacement parts. Party A(buyer) bear two-way shipping freights.

- (3) Warranty: (Please kindly note: warranty not include artificial factors breaking, such as laser generator broke by water quality).
- A. 808nm diode-300W: 1 year for mainframe, 5 millions shots warranty for handpiece.
 - B. 808nm diode-500W: 2 years for mainframe, 10 millions shots warranty for handpiece.
 - C. 808nm diode-1200W: 1 year for mainframe, 2 years with unlimited shots warranty for handpiece.
 - D. We offer technical support for lifetime for all machines.
- (4) Within 15 days from the date of machine shipment from factory, if the machine has quality problems (Non-artificial and force majeure), part B(seller) will provide replacement services. party A (seller) should ship it back and bear two-way shipping freights.
- (5) Within warranty (15th day to one year after the machine shipment), for machine quality problems (Non-artificial and force majeure).
- A. If need to send mainframe back for repairing, Party A(buyer) should bear two-way freight. The machine has to be sent back with original packaging. party B(seller) will not be responsible for the damage caused by packaging problems.
 - B. If no need to send mainframe back, Party A(buyer) can send the accessory back to repair. For some emergency and only if machine is within three months, party B(seller) can send the new accessory to change. And the old accessory must be send back to Party B(seller) within three months. Related freight should be borne by party A(buyer).
- IMPORTANT NOTICE:** Within one month, the broken parts must be returned to party B(seller). If broken parts cannot be returned in time, it is considered to be shipped as new order spare parts. The sales representative bears 40% of the cost and the after-sales department bears 60%. The party A(buyer) is considered to give up any warranty or service for this instrument. Party A(buyer) bears the two-way shipping cost
- C. For engineer on-site service, Party A(buyer) should provide with travel fee (traffic, hotel, meals).
- (6) After the machine is shipped more than one year, if machine quality problems, Party

B(seller) provides Party A(buyer) maintenance service with the cost of replacement parts. If need to send back to factory for repairing, party A (buyer) bears two-way freights.

(7) Notes:

- A. All our products are precision instruments, cannot be carried always. If any quality problems caused by the shipping, party B(seller) would provide service, but it is not the responsibility of party B(seller).
- B. Please keep the original package for the maintain and repair.
- C. Barcode is the only factory warranty certificate, party B(seller) don't provide warranty service without the barcode. The company reserves the sale described in the final interpretation and modification rights.